



Unleash Your Potential

The Eppendorf epMotion® Series:
Discover our automated liquid handling solutions



Meet the New epMotion®!

We do not want you to work like a robot, so simply let our epMotion devices do the work for you. Automation liberates you by allowing you to concentrate on tasks that require not only your hands but first and foremost your brain.

Robust Hardware


- > epMotion systems are available in three sizes, all designed with minimal footprint
- > Self checkup to reduces startup time
- > High precision dispensing tools covering a volume range from 200 nL to 1,000 µL
- > Optical sensor verifies worktable setup before the run



[Read more](#)

epBlue – Intuitive Software


- > Enjoy the drag-and-drop-based programming
- > Learn to program and adjust methods in hours, not days
- > Optional software extensions for handling barcodes and supporting regulated environments compliant to e. g. 21 CFR part 11

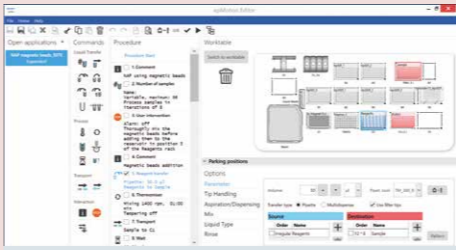

[Read more](#)

epMotion Accessories, Consumables & Services

- > Get maximum flexibility with our portfolio of epMotion accessories
- > As an open platform epMotion supports a large variety of labware
- > High quality consumables guarantee high accuracy and reproducibility of your results
- > Choose our premium service to support and maintain your equipment operation

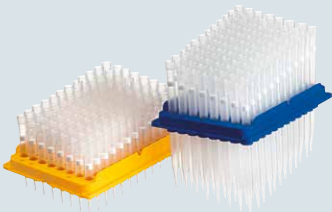

[Read more](#)

Let us focus on what really matters: Improving the world through science

This is why we get up every day, go to the lab and make new discoveries!
Let us take our time to create new experiments and analyze their results
instead of doing tedious routine pipetting work.

Works Hard, so You Can Work Smart


Before automation, researchers (like you) spent most of their working hours handling samples. That meant pipetting. A lot. Day after day.

Today is a new day. Use your valuable time for what really matters. Create experiments and evaluate results.

Replaces Your Hands, Not Your Brain

epMotion takes this (pain)work off your hands. It pipettes precisely and prevents you from developing fatigue-related injuries like Repetitive strain injury (RSI).

Never worry about sample quality again. epMotion is the friendly robot you can trust. It does its job reliably from the first to the last sample. Change the world, not your posture. Pipetting has never been this ergonomic.

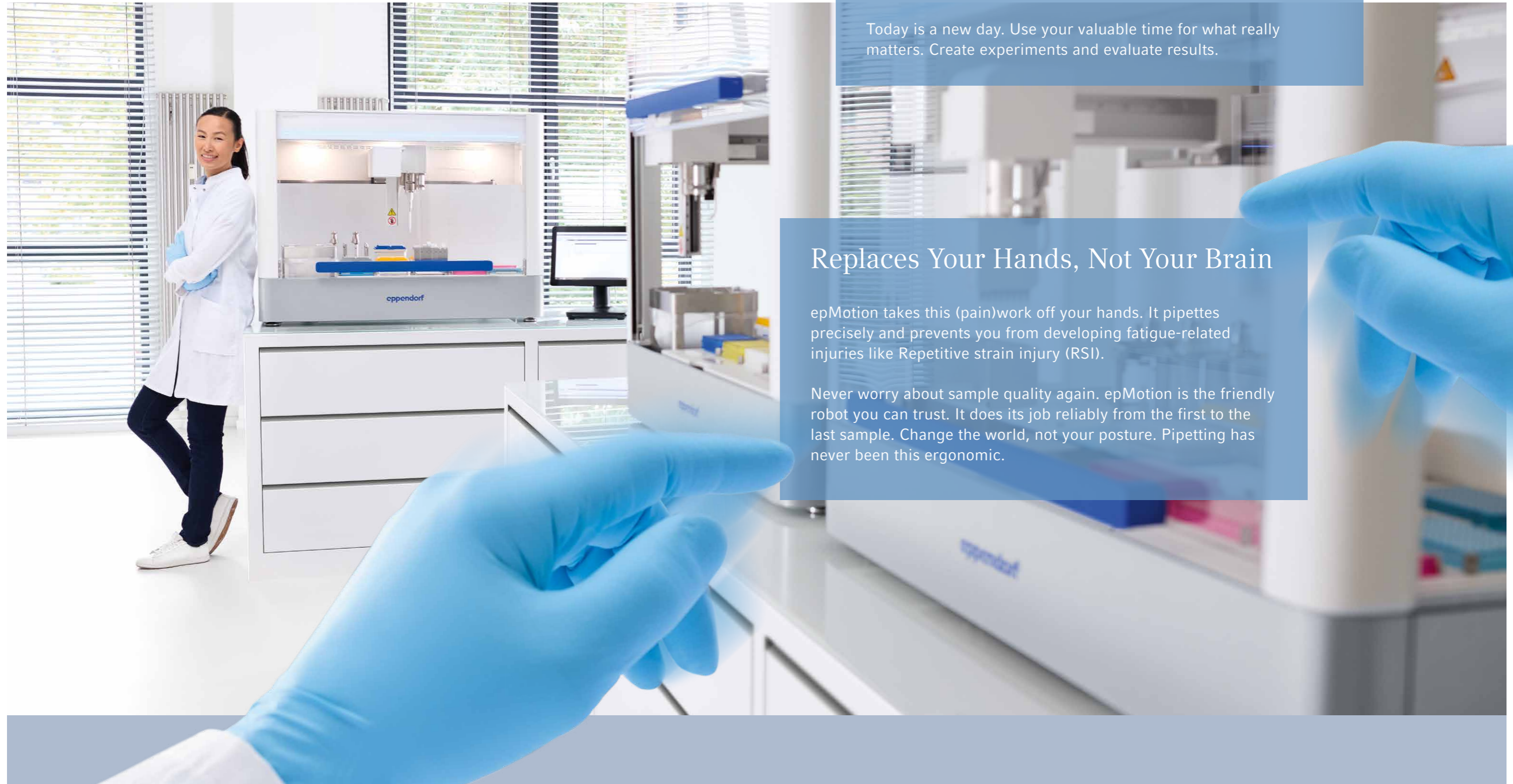

Works like a Robot, so You Don't Have to

The epMotion has its limits. It cannot think like you. But it does not forget any pipetting steps. It processes small volumes that are barely visible to the human eye with precision.

Be strong. You cannot pipette like a robot. But you are smarter. Join forces with the epMotion to take the quality of experiments to a whole new level.

Operates Tirelessly Without Getting Tired

And there is something else the epMotion does not have to do. Sleep! It never gets tired. That is why it achieves a much higher throughput.

With your research, you are expanding the boundaries of humankind. It is more than OK to get tired. But please get tired because you have been able to do what is important to you. Work with joy.

Doesn't Want Your Job, Just Your Work

Let us talk about joy. epMotion loves your tried-and-tested protocols. It even likes to track every method step and provides a detailed protocol for you.

1000

Are you a fan of painstaking notation? Then this might be disappointing for you: The epMotion does all the documentation for you during its use. Do not be sad, take a break instead!

Flexible, Yet Intuitive!

Your sample journey is different from case to case. epMotion addresses your needs, giving you full flexibility and endless combination possibilities. In just a few clicks, you decide which consumables are used and whether shaking or heating is required.

Adjust sample numbers, pipetting volumes and incubation times quickly and intuitively before the start of every protocol. No programming skills required. Just exactly what you need: you are a natural scientist not a robotics engineer.

Eppendorf epMotion® Portfolio



epMotion® 96

The Eppendorf epMotion 96 is a semi-automated, 96-well multichannel pipette for fast and precise parallel microplate processing. It is an easy to use bench top liquid handling system for high precision pipetting in 96 and 384 well plates. Without changes to the system, a large volume range of 0.5 µL to 300 µL is available for convenient use over a range of applications. Its ergonomic design and intuitive handling makes the epMotion 96 a great tool for anyone in the lab who needs fast and precise liquid handling in 96 format.

Product Feature Highlights epMotion 96/96xl

- > Electronic pipetting with parallel piston movement for better precision and reproducibility
- > Two tip sizes covering the 0.5 µL to 1,000 µL operational range for maximum accuracy
- > Auto-detection of tip sizes without the need to change heads
- > Use 96 tips at once or use tips column-wise (8, 16, 24, etc.)



epMotion® 5070

Our smallest member of the epMotion family is the most compact solution for accurate and reproducible automated pipetting. This makes the epMotion 5070 a perfect match for any routine application such as PCR and qPCR setup, serial dilutions, reagent distribution, sample transfer from tubes to plates, and sample normalization.

Product Feature Highlights epMotion 5070

- > 4-position worktable
- > Volume range 0.2 to 1,000 µL (depending on used dispensing tools)
- > Automatic tool exchange for 2 tools
- > Small footprint of 65 × 48 cm fits small lab benches
- > MultiCon PC controller with simulation, network and software upgrade options
- > Intuitive drag-and-drop-based software

Available as PCR Solution bundle with dispensing tool and accessories.



epMotion® 5073

These automated pipetting systems are perfectly suited for PCR and qPCR setup, nucleic acid purification as well as low-throughput NGS library preparation. It yet retains the flexibility that allows the use as open platforms for diverse liquid handling tasks. The epMotion 5073 series automates and simplifies what are traditionally complex, labor-intensive pipetting tasks, saving time and improving the reproducibility of results.

Product Feature Highlights epMotion 5073

- > 6-position worktable
- > Volume range 0.2 to 1,000 µL (depending on used dispensing tools)
- > Automatic tool exchange for 3 tools
- > Option for gripper, one thermal module*¹ or Eppendorf ThermoMixer®*²
- > Optional UV lamp and air filter system for decontamination and clean air condition
- > MultiCon PC controller with simulation, network and software upgrade options
- > Intuitive drag-and-drop-based software

Available as NGS Solution bundle with specific software license, dispensing tools, accessories and consumables.

^{*1} Available only with 5073i
^{**2} Already included 5073t



epMotion® 5075

The epMotion 5075 is the ideal solution for various liquid handling demands. It offers the same outstanding accuracy and precision as epMotion 5070 & 5073. The available options make the 5075 an excellent and most flexible device for applications such as NGS library preparation, PCR and qPCR setup, magnetic-bead-based and filter-based purification, cell-based assays or any routine pipetting tasks.

Product Feature Highlights epMotion 5075

- > Up to 15 worktable positions
- > Volume range 0.2 to 1,000 µL (depending on used dispensing tools)
- > Automatic tool exchange for 4 tools
- > Option for gripper, 1-3 thermal modules, Eppendorf ThermoMixer® and vacuum station
- > Optional UV lamp and air filter system for decontamination and clean air condition
- > MultiCon PC controller with simulation, network and software upgrade options
- > Intuitive drag-and-drop-based software

Available as NGS Solution bundle with specific software license, dispensing tools, accessories and consumables.



Eppendorf epMotion® 96 and 96xl

The epMotion 96/96xl is a semi-automated, 96-well multichannel pipette for fast and precise parallel microplate processing. It is an easy to use bench top liquid handling system for high precision pipetting in 96 and 384 well plates. Without changes to the system, a large volume range is available for convenient use over a range of applications. Its ergonomic design and intuitive handling makes the epMotion 96 a great tool for anyone in the lab who needs fast and precise liquid handling in 96 format.

Applications

- > Replication and reformatting of microplates
- > Cell seeding and media change
- > Reagent and compound addition
- > Magnetic-bead-based clean-ups
- > Cell-based assays
- > ELISA handling (plate coating and washing)
- > Biochemical assays

- > Volume ranges 0.5 µL to 300 µL (96) or 5 µL to 1,000 µL (96xl)
- > Electronic pipetting with synchronous 96-channel piston movement for better precision and reproducibility of results
- > Auto-detect function for tip size
- > Use 50 µL, 300 µL and 1,000 µL reload tips for best accuracy
- > Use tips column-wise (8, 16, 24, etc.) for serial dilutions and similar applications
- > Intuitive and industry proven software concept and convenient touch screen control
- > Intelligent, preset applications: aspiration, dilution, multi dispense, pipette and mix
- > Available with 1 position or 2-position slider
- > Compact design to fit laminar flow hoods
- > Reduced risk of repetitive strain injury (RSI)

Ease Your Workflow



Optional: Remove tips or add tips from the reload rack (always work with full columns of tips).

1. Insert the appropriate tip type for 0.5–50 µL or 5–300 µL or 10–1,000 µL into the loading frame.

2. Slide loading frame into epMotion 96 and apply lever to easily attach tips to the 96 nose-cone array.

3. Select mode, e.g. »Pipette + Mix«, and adjust plate on lifting table towards the tips (if 2-position slider is mounted, source or destination position is selected).

4. Touch »Aspirate« to take liquid from source plate, change to destination plate and touch »Dispense«.



Please visit our website:
www.eppendorf.com/epMotion96

epMotion® 5070



The epMotion 5070 is our most compact solution for accurate and reproducible automated pipetting. With four worktable positions and a small footprint the epMotion 5070 is a perfect match for any routine application such as PCR and qPCR setup, serial dilutions, reagent distribution, sample transfer from tubes to plates, and sample normalization.

Applications

- > PCR/qPCR setup
- > Tube to plate transfer
- > Plate reformatting
- > Serial dilutions
- > Cherry picking and Normalization
- > Sample and reagent transfers

- > 4 SLAS/ANSI deck positions
- > Volume range 0.2 to 1,000 µL (depending on used dispensing tools)
- > Automatic tool exchange for 2 tools
- > Compatible with micro- and PCR plates up to 384 wells
- > Compatible with tubes from 0.2 mL to 50 mL
- > Optical sensor for detecting liquids, labware, tips
- > Liquid level detection without conductive tips
- > Fully enclosed, compact housing
- > MultiCon PC controller with simulation, network and software upgrade options
- > Intuitive drag-and-drop-based software
- > Small footprint of 65 × 48 cm fits small lab benches
- > Available as PCR Solution bundle with dispensing tool and accessories



Please visit our website:
www.eppendorf.com/epMotion5070

epMotion® 5073I



Why spend valuable time on tedious liquid handling procedures when you can achieve fast, accurate and reproducible results with the **epMotion 5073I** automated liquid handler? Perfect for medium-throughput setup for PCR and qPCR, DNA/RNA normalization, cherry picking, dilution series and reformatting, the epMotion 5073I helps you perform pipetting tasks with ease and greater efficiency.

Applications

- > Assay setup
- > Simple cell culture applications
- > PCR/qPCR setup
- > Tube to plate transfer
- > Plate reformatting
- > Serial dilutions
- > Cherry picking and Normalization
- > Sample and reagent transfers
- > Sample mixing and tempering

- > 6 SLAS/ANSI deck positions,
- > Volume range 0.2 to 1,000 µL (depending on used dispensing tools)
- > Automatic tool exchange for 3 tools and gripper
- > Compatible with micro- and PCR plates up to 384 wells
- > Compatible with tubes from 0.2 mL to 50 mL
- > Optical sensor for detecting liquids, labware, tips
- > Liquid level detection without conductive tips
- > Fully enclosed, compact housing
- > MultiCon PC controller with simulation, network and software upgrade options
- > Intuitive drag-and-drop-based software
- > Optional UV lamp and air filter system for decontamination and clean air conditions
- > Optional thermal module for heating or cooling of samples or reagents
- > Integrated LEDs for deck illumination and visual feedback of system status like idle, run in progress or run finished
- > Optional software extensions for handling barcodes and supporting regulated environments compliant to e. g. 21 CFR part 11



Please visit our website:
www.eppendorf.com/epMotion5073I

epMotion® 5073t



The epMotion 5073t automated liquid handling system allows you to perform fast and efficient NGS library preparation of 1 to 24 samples. It is also ideally suited for automating more complex, time consuming liquid handling tasks such as DNA/RNA extraction. Save time and improve the reproducibility of your results with the epMotion 5073t!

Applications

- > DNA and RNA purification
- > NGS library preparation (low sample numbers)
- > PCR/qPCR setup
- > Tube to plate transfer
- > Plate reformatting
- > Serial dilutions
- > Cherry picking and Normalization
- > Sample and reagent transfers
- > Sample mixing and tempering

- > 6 SLAS/ANSI deck positions
- > Volume range 0.2 to 1,000 µL (depending on used dispensing tools)
- > Automatic tool exchange for 3 tools and gripper
- > Compatible with micro- and PCR plates up to 384 wells
- > Compatible with tubes from 0.2 mL to 50 mL
- > Optical sensor for detecting liquids, labware, tips
- > Liquid level detection without conductive tips
- > Fully enclosed, compact housing
- > MultiCon PC controller with simulation, network and software upgrade options
- > Intuitive drag-and-drop-based software
- > Optional UV lamp and air filter system for decontamination and clean air conditions
- > Integrated Eppendorf ThermoMixer with ^{2D}Mix-Control technology
- > Integrated LEDs for deck illumination and visual feedback of system status like idle, run in progress or run finished
- > Optional software extensions for handling barcodes and supporting regulated environments compliant to e. g. 21 CFR part 11
- > Available as NGS Solution bundle with specific software license, dispensing tools, accessories and consumables



Please visit our website:
www.eppendorf.com/epMotion5073t

epMotion® 5075I



Do not stress yourself by doing complex pipetting tasks when you can perform fast, accurate and reproducible liquid handling operations with the **epMotion 5075I** automated liquid handler. The epMotion 5075I offers 15 worktable positions and is well suited for high-throughput preparation of plates for PCR and qPCR, DNA/RNA normalization, cherry picking, dilution series, reformatting and immunoassays. Say goodbye to tedious pipetting steps!

Applications

- > Assay setup
- > Simple cell culture applications
- > PCR/qPCR setup
- > Tube to plate transfer
- > Plate reformatting
- > Serial dilutions
- > Cherry picking and Normalization
- > Sample and reagent transfers
- > Sample mixing and tempering

- > 15 SLAS/ANSI deck positions
- > Volume range 0.2 to 1,000 µL (depending on used dispensing tools)
- > Automatic tool exchange for 4 tools and gripper
- > Compatible with micro- and PCR plates up to 384 wells
- > Compatible with tubes from 0.2 mL to 50 mL
- > Optical sensor for detecting liquids, labware, tips
- > Liquid level detection without conductive tips
- > Fully enclosed, compact housing
- > MultiCon PC controller with simulation, network and software upgrade options
- > Intuitive drag-and-drop-based software
- > Optional UV lamp and air filter system for decontamination and clean air conditions
- > Optional up to 3 thermal modules for heating or cooling of samples or reagents
- > Integrated LEDs for deck illumination and visual feedback of system status like idle, run in progress or run finished
- > Optional software extensions for handling barcodes and supporting regulated environments compliant to e. g. 21 CFR part 11



Please visit our website:
www.eppendorf.com/epMotion5075I

epMotion® 5075t



The **epMotion 5075t** automated liquid handling system allows you to perform high-throughput NGS library preparation using magnetic beads with speed and efficiency. It is also ideally suited for automating more complex, time consuming liquid handling tasks such as DNA/RNA extraction and ELISA. Save time and improve the reproducibility of your results with the epMotion 5075t!

Applications

- > DNA and RNA purification
- > NGS library preparation (medium to large sample numbers)
- > PCR/qPCR setup
- > Tube to plate transfer
- > Plate reformatting
- > Serial dilutions
- > Cherry picking and Normalization
- > Sample and reagent transfers
- > Sample mixing and tempering

- > 14.5 SLAS/ANSI deck positions
- > Volume range 0.2 to 1,000 µL (depending on used dispensing tools)
- > Automatic tool exchange for 4 tools and gripper
- > Compatible with micro- and PCR plates up to 384 wells
- > Compatible with tubes from 0.2 mL to 50 mL
- > Optical sensor for detecting liquids, labware, tips
- > Liquid level detection without conductive tips
- > Fully enclosed, compact housing
- > MultiCon PC controller with simulation, network and software upgrade options
- > Intuitive drag-and-drop-based software
- > Optional UV lamp and air filter system for decontamination and clean air conditions
- > Integrated Eppendorf ThermoMixer with ^{2D}Mix-Control technology
- > Optional up to 2 thermal modules for heating or cooling of samples or reagents
- > Integrated LEDs for deck illumination and visual feedback of system status like idle, run in progress or run finished
- > Optional software extensions for handling barcodes and supporting regulated environments compliant to e. g. 21 CFR part 11
- > Available as NGS Solution bundle with specific software license, dispensing tools, accessories and consumables



Please visit our website:
www.eppendorf.com/epMotion5075t

epMotion® 5075v



Complex applications, such as nucleic acid purification or solid phase extraction, should not tie you to the workbench. With integrated vacuum filtration and a combined mixing and temperature module (only for 5075vt) in one system, the **epMotion 5075v** and **epMotion 5075vt** allow true walk-away automation of complex purification or extraction protocols. Increase your lab productivity and free up time for other tasks!

Applications

- > DNA and RNA purification
- > Solid phase extraction
- > PCR/qPCR setup
- > Tube to plate transfer
- > Plate reformatting
- > Serial dilutions
- > Cherry picking and Normalization
- > Sample and reagent transfers
- > Sample mixing and tempering

- > 12 SLAS/ANSI deck positions
- > Volume range 0.2 to 1,000 µL (depending on used dispensing tools)
- > Automatic tool exchange for 4 tools and gripper
- > Compatible with micro- and PCR plates up to 384 wells
- > Compatible with tubes from 0.2 mL to 50 mL
- > Optical sensor for detecting liquids, labware, tips
- > Liquid level detection without conductive tips
- > Fully enclosed, compact housing
- > MultiCon PC controller with simulation, network and software upgrade options
- > Intuitive drag-and-drop-based software
- > Integrated vacuum pump – silent operation, no tubing, wiring, and reservoirs to maintain
- > Integrated Eppendorf ThermoMixer with ^{2D}Mix-Control technology (only 5075vt)
- > Optional UV lamp and air filter system for decontamination and clean air conditions
- > Optional up to three thermal modules for heating or cooling of samples or reagents
- > Integrated LEDs for deck illumination and visual feedback of system status like idle, run in progress or run finished
- > Optional software extensions for handling barcodes and supporting regulated environments compliant to e. g. 21 CFR part 11



Please visit our website:
www.eppendorf.com/epMotion5075vt



epMotion® 5070/5073/5075

Feature Overview




Optical Sensor
Checks worktable loading for tips (type and number), labware type, and even liquid levels in vessels for safe operation.



Solid and Liquid Waste System
Extended waste capacity for solid and liquid waste. Easy and safe handling and emptying of both types of waste.



Clean Cap
with UV light for decontamination and air filter system for clean air inside. Optional for epMotion 5073 and 5075 versions.

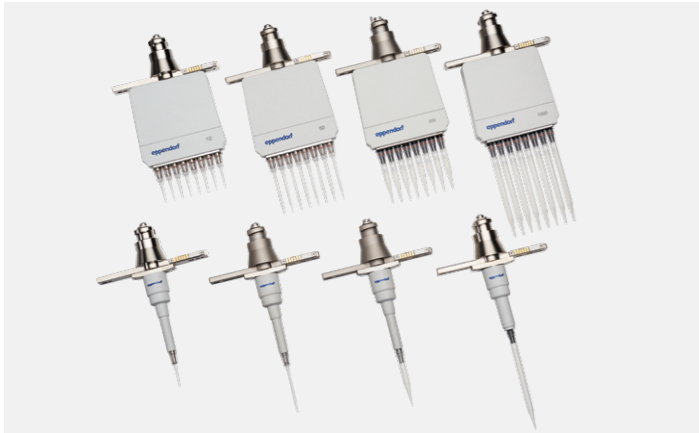


Dispensing Tools
8 autoclavable dispensing tools from 200 nL to 1,000 µL are available for programming; single- or multi-channel; all calibrated for utmost precision.



Gripper Tool
Allows transport and stacking of labware on worktable of 5073 or 5075 versions.

Accurate Pipetting Performance
Achieve highest pipetting precision and accuracy with free-jet or contact dispensing using the Eppendorf air-cushion pipetting technology together with dedicated epT.I.P.S. Motion in a volume range from 0.2 to 1,000 µL. While free-jet pipetting minimizes cross contamination, contact pipetting can be used for better precision of very small volumes.

Use Your Favorite Labware
Take advantage of our growing database of over 1,200 labware files, which supports the labware you are already using. Our database includes a range of single tubes, from 0.2 mL PCR tubes to 50 mL conical tubes, as well as plates from 6 to 384 well format. Additionally, we offer a selection of adapters, including cooling and heating options, for microplates and deepwell plates from various suppliers.

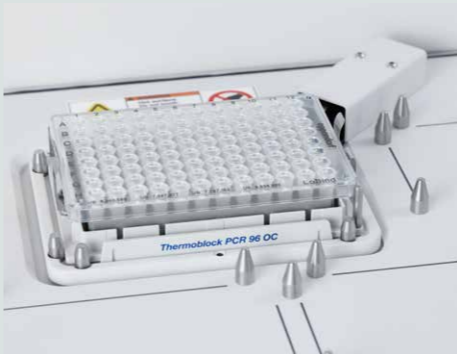

Eppendorf ThermoMixer®
Fully integrated Eppendorf ThermoMixer with 2D Mix-Control improves yield and application performance, e.g. for resuspending bacteria pellets or magnetic beads. Peltier element allows heating and cooling of sample and/or reagents.




Integrated Vacuum Station
Loaded and unloaded with the gripper, the vacuum station adapts automatically to different filter plate types. Silent operation of internal pump with no need for extra tubes, wiring or reservoir maintenance.



Status Display and Deck Illumination
LED status light bar for visual feedback of the system status and separate LED light for worktable illumination for epMotion 5073 and 5075.



Automatic Tool Exchange
Choose the tools you need for your application from 8 dispensing tools and a gripper. Dependent on the epMotion version, 2 to 4 dispensing tools can be provided on the worktable for automatic exchange.

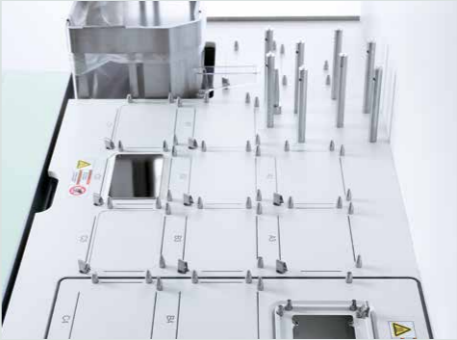

Small Footprint – Large Deck
With 4, 6, 12 or up to 15 positions in SLAS/ANSI format the epMotion family gives you all the flexibility you need for your application.

Easy & Fast Access
Control epMotion by touch or mouse for quick and easy method setup. Fast initialization of the instrument let you start right away without manual teaching or calibration of robotic system.

Housing
Completely contained housing including door safety mechanism.

Thermal module
Option for heating or cooling of samples and reagents (0–110 °C)

NGS Made Easy

NGS sample preparation is a labor-intensive process. With its precision pipetting capabilities, the epMotion reduces the hands-on time and generates high-quality NGS libraries. Eppendorf offers pre-optimized and manufacturer-qualified NGS reagent kit methods, allowing for minimal user

intervention and setup time, even for runs with low sample numbers. The system's ease of use is further enhanced by Eppendorf's expert applications support and service, ensuring successful automation of a variety of kits.



For consistent high-quality results, each workflow requires careful optimization, customization, and thorough training. To make the process even smoother, Eppendorf provides in-lab applications implementation services, making it quick and simple to add new methods to the epMotion. With this level of support, the epMotion eliminates the risk of pipetting errors, provides reproducible results, and increases overall productivity.


In addition, the epMotion NGS solutions offer users the choice for a complete solution for processing up to 24 or up to 96 samples without user-intervention, freeing up time for other tasks. The sequencing results obtained through epMotion automation are comparable to or even better than those obtained through manual preparation, making the epMotion a trustworthy solution for automating NGS library preparation.

epMotion® Software – Simplify Your Programming

Versatile Software Gives You Control

The epBlue software features an intuitive interface that allows for fast mastery and rapid method design while maintaining versatility and usability. The easy to use drag-and-drop visual interface makes setting up new methods fast and simple. A variety of smart commands are available each with its own set of preset parameters that the user can modify to fit particular needs. Comments give the option to leave notes when multiple users work on the same protocol. Features such as the ability to reuse tips and a multidispense mode minimize tip consumption leading to reduced costs per sample can be activated with a few clicks. The optional barcode reading tool allows for easy tracking of samples throughout your whole process.

- > Intuitive drag-and-drop programming for easy method setup
- > Copy/Cut and paste feature allows for fast method modification and adaptation
- > 3D view of workstation and run simulation for easier method optimization
- > Create, edit and simulate methods on any Windows 10 PC with epMotion Editor software
- > User management with different user levels enables structured use in a shared environment
- > The ability to leave comments makes it easy to collaborate on methods
- > Surface teaching allows for the user to identify the bottom tolerance of labware
- > Unlimited virtual deck positions increases the versatility of the system
- > Add Enhanced Feature Set 1 and 2 for even more features

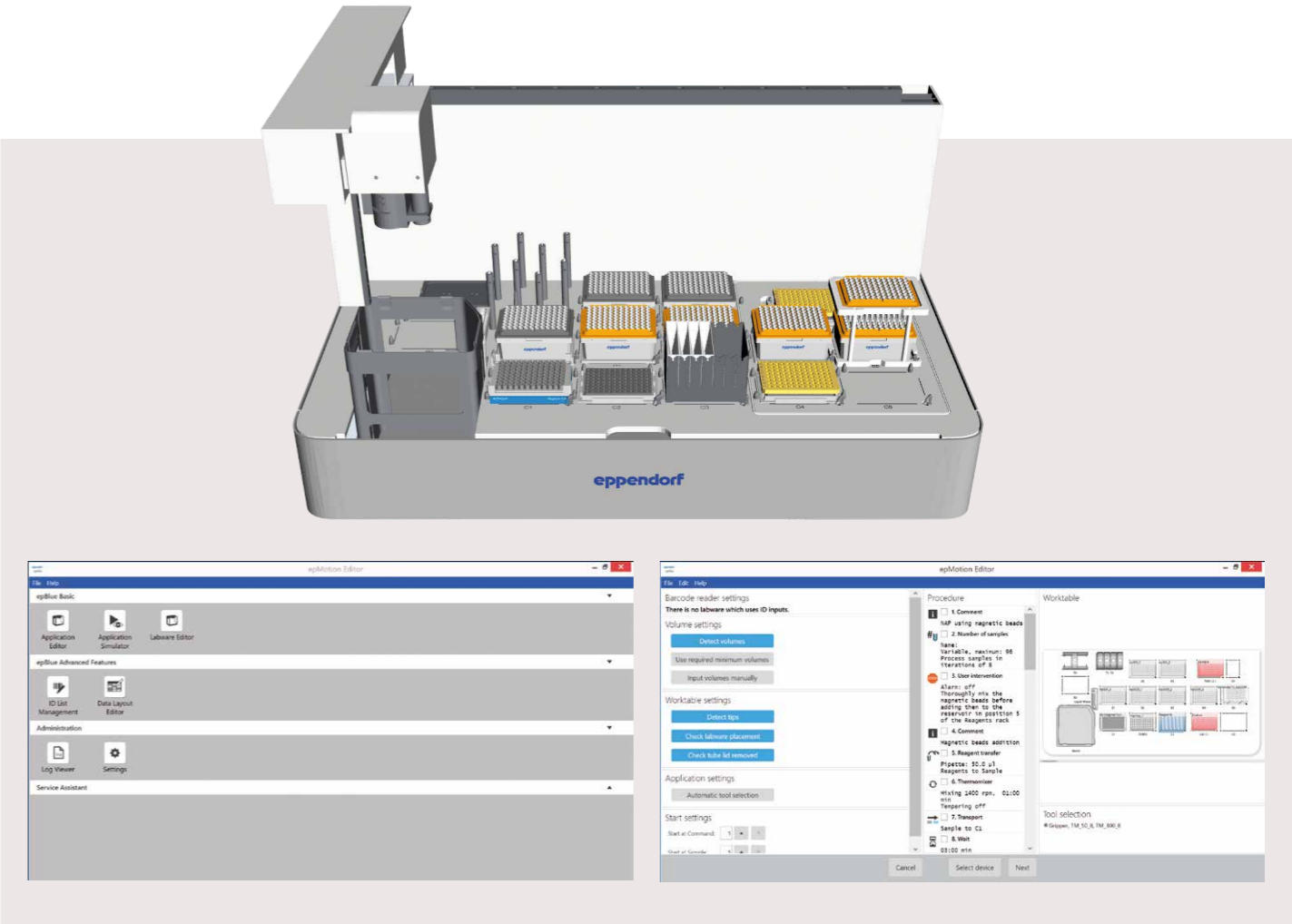

Please visit our website:
www.eppendorf.com/epmotionmethods/


Upgrades for the Features You Need

epBlue GxP – Software

The epBlue GxP software was developed to support organizations in validated regulated environments. The epBlue GxP software is designed to the properties required for processes compliant to the requirements of e.g. 21 CFR part 11 containing the following functionalities:

- > Electronic records and signatures
- > User level management and access control
- > Audit trail and log file
- > Revision management
- > Export and archiving of digitally signed documents
- > epBlue ID tracking using bar codes (optional extension)


In addition to the epBlue GxP software, GxP IQ and GxP OQ services are available for all epMotion models considering the currently valid regulations and standards to meet customer needs in terms of compliance.



epBlue ID – Secure barcode scanning-tracking-documentation

epBlue ID software module allows for safe data exchange with laboratory information management systems (LIMS) and simplifies external communication. The barcode scanner enables the user to record barcodes on all tube types and microplates. Reagents can be documented with type and lot number. After manually scanning the barcodes, epBlue ID will store the IDs in its database from which data can be verified at any time. A result file containing sample IDs and their final location is generated by epBlue ID when the liquid handling process is completed. The result ID list can then be exported to a network drive and uploaded to a LIMS.

- > Visual guidance of the scanning process
- > Documentation of reagent type and batch
- > Compatible with LIMS generated worklists

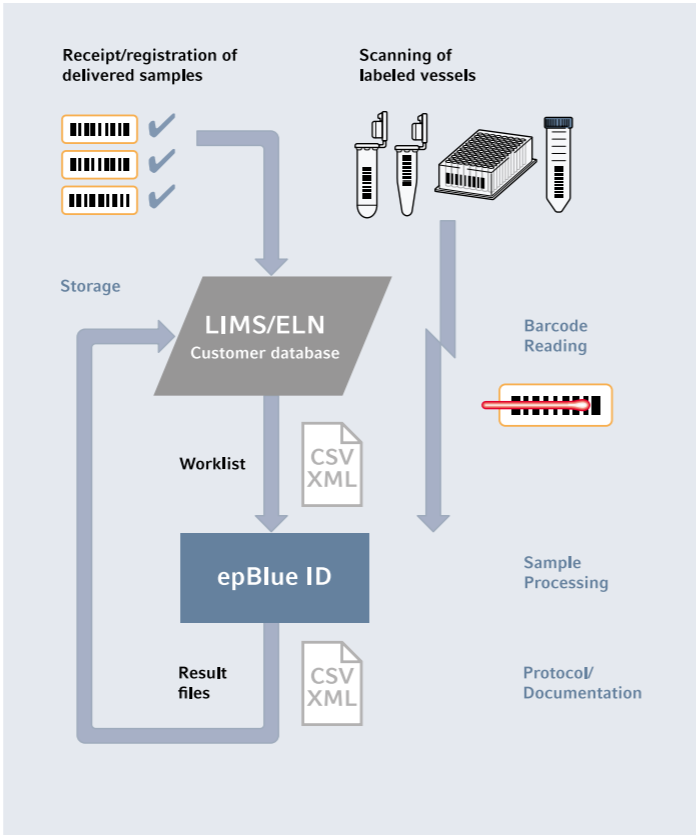


Please visit our website:
www.eppendorf.com/epBlue

Control at Your Fingertips

The Eppendorf MultiCon PC versions define a new dimension in ease-of-use and safety. As premium product, the industrial PC ensures maximum reliability by standardization and has undergone multiple quality checks before delivery. The MultiCon controller fulfills the high expectations for best performance for the epMotion.



- > Large touch screen and mouse and keyboard control
- > Network interface allows the remote control of epMotion via direct Ethernet connection
- > Pre-installed epBlue software
- > CSV and XML file exchange via USB or LAN
- > Safe operation due to data-base repository
- > 5 USB & 1 Ethernet ports
- > Optional software extensions for handling barcodes and supporting regulated environments compliant to e.g. 21 CFR part 11



Accuracy and Reproducibility Delivered

epMotion dispensing tools

- Based on Eppendorf air cushion pipetting technology
- > Free-jet pipetting, no tubes, no wires, no air bubbles to minimize cross-contamination
 - > Use contact pipetting for small volumes, from 1 µL down to 200 nL
 - > Excellent pipetting precision, below 2 % CV at 1 µL*
 - > Volume range: 0.2–10 µL; 1–50 µL; 20–300 µL and 40–1,000 µL
 - > Single or 8-channel
 - > Autoclavable & easy to maintain
 - > Automatic calibration warning

* Eppendorf application note 168:
»Evaluation of the Eppendorf epMotion® Pipetting Tools using the Artel MVS®.«



epT.I.P.S.® Motion

- Special tips for use in automated systems. Each tip is inspected for straightness prior to packaging. This ensures a perfect handling of 96- and 384-well micro- and PCR plates
- > Color-coded tray for easy volume identification
 - > Non-carbonized, pure virgin polypropylene, both tip and box are recyclable
 - > Standard, filter, and sterile tips available
 - > PCR clean quality grade for filter tips
 - > Also available as eco-friendly reload system (ACT certified)



Re-use tips feature

- The new epMotion systems can re-use tips as e.g. economical choice for multiple wash steps of the same well. Here special epT.I.P.S. Motion SafeRack tips are recommended
- > Individual compartments to avoid cross-contamination between used tips
 - > Available as standard and PCR clean quality



Services for epMotion® – Keep Your Process in Motion

It’s all about accuracy and reliability

The epMotion and its dispensing tools are precision instruments. As such, they require regular maintenance, which help prevent failures and keep your processes moving forward. Regular calibration of the dispensing tools will ensure system accuracy, precision and reliability.

Our Service and Extended Warranty Agreement options will allow you to maintain premium functioning of your epMotion over the full lifetime of the system. Eppendorf’s highly skilled service technicians perform the adjustments and repairs to return your system to manufacturer specifications.

Services include:

- > Cleaning, inspection and system maintenance
- > Eppendorf quality spare parts
- > Software update as needed
- > Instrument adjustments as needed
- > Functional tests with service diagnostic methods
- > Calibration of dispensing tools
- > Calibration of level sensor, carrier arm, thermal module and vacuum manifold
- > Operational Qualification certification
- > Full documentation
- > 1 year extended warranty

epMotion®	96	5070	5073	5075
Installation	0082 030 171	0082 030 081	0082 030 111	0082 030 141
PremiumCare	0082 030 076	0082 030 107	0082 030 108	0082 030 075
AdvancedCare	0082 030 074	0082 030 105	0082 030 106	0082 030 073
Extended Warranty	0082 030 070	0082 030 100	0082 030 101	0082 030 069
Extended Warranty plus	0082 030 072	0082 030 102	0082 030 103	0082 030 071
Preventive Maintenance	0082 030 002	0082 030 010	0082 030 011	0082 030 012
IQ/OQ GxP	0082 030 036	0082 030 030	0082 030 031	0082 030 032
OQ GxP	0082 030 016	0082 030 040	0082 030 041	0082 030 042
Intro user training	0082 030 185	0082 030 095	0082 030 125	0082 030 155

Application Implementation		Calibrations	
Level 1	0082 030 160	epMotion® 1-ch. Tool Calib. (10 M/test vol)	0050 600 255
Level 2	0082 030 161	epMotion® 8-ch. Tool Calib. (10 M/test vol)	0050 600 259
Level 3	0082 030 162	epMotion® 1-ch. Tool Calib. (4 M/test vol)	0082 030 077
Level 4	0082 030 163	epMotion® 8-ch. Tool Calib. (4 M/test vol)	0082 030 078
		epMotion® 96 Calibration	0082 030 020

For more detailed information on our service portfolio please download the epMotion Services brochure here.




Further information:
www.eppendorf.com/epServices
For local offers please visit our local websites.
Services may vary according to country.
Contact your local Eppendorf Service organization.




All Your Favourite Labware – epMotion® Accessories


Whether 0.2 mL PCR tubes, 1.5, 2 or 5 mL test tubes, 15 or 50 mL conical tubes, six-well plates, 96/384-well microplates, PCR plates, or deep-well plates – with epMotion accessories, you can handle it.




ReservoirRack
Enables positioning of up to 7 reagent rack modules or reservoirs with with 10 mL, 30 mL and 100 mL filling volume.

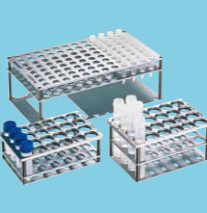


Racks for single test tubes
For micro test tubes, glass, or plastic tubes. Fifteen different formats available.




ReservoirRack 3
Small rack for epMotion 5075t for 3 reagent rack modules or reservoirs with 10 mL, 30 mL or 100 mL filling volume.




Thermoracks for 24 × Safelock 0.5/1.5/2 mL tubes
Temperature-controlled when used with a thermal module. Available also for Cryotubes and as version for Eppendorf ThermoMixer®.




ReservoirRack modules TC
10 module sizes are available to use tubes from 0.2 mL PCR, 5 mL to 50 mL tubes. Temperature-controlled when used with a thermal module.

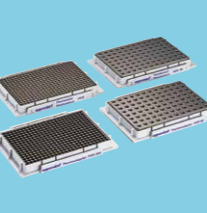


Thermoracks and adapters
For use with 96-well or 384-well plates and single 0.2 mL tubes. Thermoracks can be transported with gripper. On Thermoadapters plates can be exchanged by the gripper.




ReservoirRack Module TC 5.0 mL
For 4 × Eppendorf Tubes® 5.0 mL. Up to 3 modules can be placed in ReservoirRack, temperature control possible.




Thermoadapter DWP 96
For heating or cooling of 1 mL deep-well plates. Plates can be exchanged by the gripper.

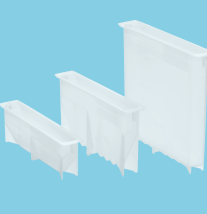


10 mL, 30 mL and 100 mL reagent reservoirs
The reagent reservoirs are PCR clean, autoclavable, and can be placed in the ReservoirRack. Both can be temperature controlled with special ReservoirRack module.

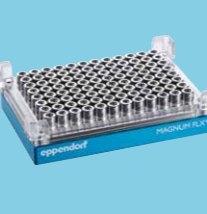


Eppendorf Magnum FLX® Magnet Adapter
universal ring magnet plate for rapid magnetic bead separation. Supports most 96-well plates (standard, PCR, deepwell).




Height adapter
Three height adapters enable exact height level adjustment and accelerated processing of microplates.




Vac Thermo lid
Use with Vacuum station to apply heated air to samples for improved drying.




Gripper
For transport of plates, thermoblocks and other labware on worktable of epMotion 5073 & 5075.

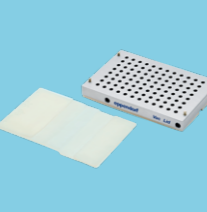


Vacuum lid
For use with Vacuum station to apply vacuum to a reduced number of samples. Grippable and special sealing mats for lid available.




TipHolder 73
enables positioning epT.I.P.S.Motion Reload trays above a rack of epT.I.P.S. Motion and thus doubles the tip capacity on a deck position.

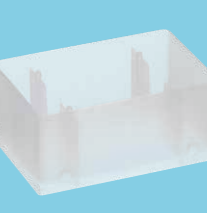


Reservoir 400 mL
Large-volume container in SLAS/ANSI format as e.g. liquid waste collection vessel for vacuum applications on the epMotion, autoclavable.




TipHolder
Autoclavable, aluminum Tipholder for tip tray of reload system, reduces waste by up to 40 percent.



Liquid Waste System
System for easy and safe handling of liquid waste.



Reservoir Rack Module NGS
for use with epMotion ReservoirRacks, for distributing samples and reagents for NGS-specific applications, space for 20 × ILMN tubes, 4 × Eppendorf Tubes® 5.0, 4 × Eppendorf Tubes 1.5 mL.



Reservoir Rack Module PCR
for use with epMotion ReservoirRacks, for distributing samples and reagents for PCR, space for 48 × PCR tubes 0.2 mL, 12 × Eppendorf Tubes 1.5/2.0 mL and 1 × tube 5.0 mL.

Technical Specifications

Dimensions of device		5070	5073	5075
	Width*1	25.6 in/65 cm	26.8 in/68.0 cm	42.5 in/108.0 cm
	Depth	18.9 in/48 cm	24.5 in/62.3 cm	24.5 in/62.3 cm
	Height	24.8 in/63 cm	32.0 in/81.2 cm	32.0 in/81.2 cm

*1 Total width needed on bench with MultiCon and Keyboard add 19.6 in/50 cm

Weight (without packaging and accessories)					
5070	45 kg	99.2 lb	5075l	122 kg	269.0 lb
5073l	86 kg	189.6 lb	5075t	124 kg	273.4 lb
5073t	93 kg	205.0 lb	5075v	123 kg	271.2 lb
CleanCap (5073)	4.5 kg	9.9 lb	5075vt	128 kg	282.2 lb
Thermal Module (5073, 5075)	3.1 kg	6.7 lb	CleanCap (5075)	6.3 kg	13.9 lb

Power supply	5070	5073	5075
Voltage	100–240 V ±10 %	100–240 V ±10 %	100–240 V ±10 %
Frequency	50–60 Hz ±5 %	50–60 Hz ±5 %	50–60 Hz ±5 %
Standby	50 W	50 W	50 W
Max output	200 W	600 W	700 W

Pipetting performance specifications measured based on ISO 8655-6			
Single-channel tool	Volume	Systematic measurement error	Random measurement error
TS 10	10 µL	± 1.2 %	± 0.25 %
	1 µL	± 5.0 %	± 3.0 %
	200 nL*2	± 25 %	± 19.8 %
TS 50	50 µL	± 1.0 %	± 0.4 %
	1 µL	± 15 %	± 5 %
TS 300	300 µL	± 0.6 %	± 0.3 %
	20 µL	± 4 %	± 2.5 %
TS 1000	1,000 µL	± 0.7 %	± 0.15 %
	40 µL	± 5 %	± 1.5 %

in pipetting mode, free jet
*2 in pipetting mode, contact, typical pipetting performance

Typical pipetting performance determined using the Artel® MVS® (Multichannel verification system)*3			
Single-channel tool	Volume	systematic measurement error	random measurement error
TS 10*4	10 µL	± 0.37 %	± 0.27 %
	1 µL	± 1.13 %	± 0.64 %
	200 nL	± 8.16 %	± 1.76 %
TS 50	50 µL	± 0.64 %	± 0.20 %
	1 µL	± 5.43 %	± 1.58 %
TS 300	300 µL	± 0.70 %	± 0.13 %
	20 µL	± 0.85 %	± 0.20 %
TS 1000	200 µL	± 1.35 %	± 0.10 %
	40 µL	± 0.82 %	± 0.23 %

*3 Eppendorf Application Note 168: »Evaluation of the Eppendorf epMotion® Pipetting Tools using the Artel MVS®«
*4 Eppendorf Application Note 399: »Evaluation of the Eppendorf epMotion® 10 µL Pipetting Tools Using the Artel® MVS®«

Conductor	Systematic measurement error	Random measurement error
X, Y, Z Positioning:	± 0.3 mm	± 0.1 mm
Detector		
Optical confocal infrared detector	contact-free detection of liquid levels, tools used, labware worktables, tip types and quantities	
Optical sensor	liquid surface must be 90 ±3° to the vertical plane of the optical sensor	
Gripper		
Carrying capacity	1,200 g	
Vacuum unit		
Max output	35 NI/min	
Suction range	0.1–85 kPa	
Suction time	1 sec to 36 min	
Eppendorf ThermoMixer®		
Mixing frequency	Off, 300–2,000 rpm, optimized for labware and application	
Mixing amplitude	3 mm	
Temperature range	–15 °C below RT to 95 °C measured in labware at 75 % filling height	
Temperature homogeneity	≤ 10 % on the block	
Heating rate	5 °C/min, on the block	
Cooling rate	3 °C/min, on the block, above RT	
Setting range time	5 sec–120 min	
Optional thermal module		
Setting range time	1 min–120 min	
Setting range temperature	0 °C–110 °C (32°F–230°F)	
CleanCap		
Air filter	70–80 m³/h (5073); 70–120 m³/h (5075); class E10	
UV light	254 nm (UV-C)	
UV lamp power 1m distance	39 µW/cm²	

MultiCon (15“ Display)	
Processor	Intel® Core™ e i5-8500T Processor (with Intel vPro®)
Graphics	Intel UHD Graphics 630
Memory (RAM)	8 GB 2400 MHz DDR4 SO-DIMM
Storage	128 GB M.2 SATA SSD
Chipset	Intel Q370 Express
Display	15.6” diagonal, active matrix TFT LED LCD
USB Ports	> 3 × USB 3.0 (1 on display head, 2 on I/O) > 2 × USB 2.0 (on I/O)
Serial Ports	2 × RJ45
Ethernet Port	1 × LAN RJ45 (Gigabit)
Wireless	802.11 a/b/g/n/ac
Operating system	Windows® 10 IoT Enterprise LTSC 64-bit
Power Supply	> AC input voltage: 100–240 VAC > Input frequency: 50–60 Hz > Maximum output power: 180 W
Power consumption	76 W
Dimensions	386 mm × 369 mm × 221 mm
Weight	4.9 kg (unpacked)
Accessories:	NA Power Cord, EU Power Cord, Power Adapter – 180 W, Cable RJ-45 to DB-9, Cable cover locking screw, QIG

Ordering Information

Description	International Ordering no.	North America Ordering no.
epMotion® 96		
epMotion® 96 , semi-automated electronic pipette for parallel 96 channel microplate processing, 0.5–300 µL	5069 000 012	5069 000 004
epMotion® 96, with 2-position slider , semi-automated electronic pipette for parallel 96 channel microplate processing, 0.5–300 µL	5069 000 110	5069 000 101
epMotion® 96xl , semi-automated electronic pipette for parallel 96 channel microplate processing, 5–1,000 µL	5069 000 217	5069 000 209
epMotion® 96xl, with 2-position slider , semi-automated electronic pipette for parallel 96 channel microplate processing, 5–1,000 µL	5069 000 314	5069 000 306
epMotion® 5070		
epMotion® 5070 MultiCon , completely contained housing, system incl. Eppendorf MultiCon, epBlue software, keyboard, mouse, waste box, 100–240 V ±10 %/50–60 Hz ±5 %	5070 000 282	5070 007 011
epMotion® 5070 MultiCon PCR Solution , includes MultiCon PC, dispensing tool (TS 50), PCR specific accessories, 100–240 V ±10 %/50–60 Hz ±5 %	5070 000 948	
epMotion® 5073		
epMotion® 5073l , MultiCon PC, completely contained housing, epBlue software, keyboard, mouse, system for solid and liquid waste, 100–240 V/50–60 Hz	5073 000 110	
epMotion® 5073t , MultiCon PC, completely contained housing, ThermoMixer, epBlue software, keyboard, mouse, system for solid and liquid waste, 100–240 V/50–60 Hz	5073 000 111	
epMotion® 5073t NGS Solution , MultiCon PC, ThermoMixer, 3 dispensing tools, gripper, NGS specific accessories and consumables, epBlue software with Enhanced Feature Set 1, system for solid and liquid waste, 100–240 V/50–60 Hz	5073 000 112	
epMotion® 5075		
epMotion® 5075l , MultiCon PC, completely contained housing, epBlue software, keyboard, mouse, system for solid and liquid waste, 100–240 V/50–60 Hz	5075 000 041	
epMotion® 5075t , MultiCon PC, completely contained housing, ThermoMixer, epBlue software, keyboard, mouse, system for solid and liquid waste, 100–240 V/50–60 Hz	5075 000 042	
epMotion® 5075v , MultiCon PC, completely contained housing, vacuum system with accessories, gripper, epBlue software, keyboard, mouse, system for solid and liquid waste, 100–240 V/50–60 Hz	5075 000 043	
epMotion® 5075vt , MultiCon PC, completely contained housing, vacuum system with accessories, gripper, ThermoMixer, epBlue software, keyboard, mouse, system for solid and liquid waste, 100–240 V/50–60 Hz	5075 000 044	
epMotion® 5075t NGS Solution , MultiCon PC, ThermoMixer, thermal module, 4 dispensing tools, gripper, NGS specific accessories and consumables, epBlue software with Enhanced Feature Set 1, system for solid and liquid waste, 100–240 V/50–60 Hz	5075 000 045	
Power supply cables		
EU Power Cord , with EU-plug and C13 coupling, 2.5 m	0013 563 934	
US Power Cord , with US-plug and C13 coupling, 2.5 m	0013 563 942	022377183
UK Power Cord , with UK-plug and C13 coupling, 2.5 m	0013 594 490	
China Power Cord , with China-plug and C13 coupling, 2.5 m	0013 613 952	
Australia Power Cord , with Australia-plug and C13 coupling, 2.5 m	0013 592 454	
Argentina Power Cord , with Argentina-plug and C13 coupling, 2.5 m	0013 613 973	
Software and upgrade options		
Enhanced Feature Set 1 , license for epBlue features: flexible number of samples, intelligent selection 1- or 8-channel tool, start methods at any step, e-mail notifications, requires service visit and epBlue >40.6	5075 000 964	
Enhanced Feature Set 2, license for epBlue feature: Normalization, requires service visit and epBlue 40.8 or higher	5075 000 981	
epBlue ID software and hardware upgrade set , for MultiCon versions (epMotion® SN > 4,000), barcode support include software, barcode reader and stand, not compatible with epMotion® panel or EasyCon versions	5075 002 701	
Thermal Module Upgrade Set , retrofitting of one thermal module for heating or cooling (0–110 °C) of plates, thermoadapters, thermoblocks and thermoracks	5075 757 508	
Device Hardware Upgrades and Software , Upgrade set to MultiCon, changes existing epMotion® EasyCon to PC controller. Includes PC, keyboard, mouse, epBlue SW, PCR & Prep assistant, cables	5075 000 725	
CleanCap Upgrade Set for epMotion® 5073 , includes UV lamp and air filter system for decontamination and pipetting under clean air conditions, for upgrading epMotion® 5073 with SN > 7000 and epBlue version > 40.7	5073 001 333	
CleanCap Upgrade Set for epMotion® 5075 , includes UV lamp and air filter system for decontamination and pipetting under clean air conditions, for upgrading epMotion® 5075 with SN > 10,000 and epBlue version > 50.0	5075 001 888	
Loading frame for epMotion® 96	5069 072 005	
Adapter 384-well plate for epMotion® 96	5069 073 001	
Modification kit for 2-position slider (4.5 mm increments) including constant-force spring, shielding plates (epMotion® 96 and 96xl), screws, guiding aid and holding devices (epMotion® 96 as of serial number 2008, all epMotion 96xl)	5069 878 506	
O-ring tool set , with separator, blade, hook, mounting cone, mounting cylinder	5069 858 505	
O-ring black , 100 pieces	5069 867 504	

Description	International Ordering no.	North America Ordering no.
Dispensing tools		
Highly precise pipetting heads, for use in the tool holder of the epMotion® workstation. Each dispensing tool is completely autoclavable at 121 °C, 1 bar for 20 min. A quality certificate for the measurement results accompanies each tool.		
TS 10, single-channel dispensing tool for the volume range 0.2–10 µL	5280 000 100	
TS 50, single-channel dispensing tool for the volume range 1–50 µL	5280 000 010	960001010
TS 300, single-channel dispensing tool for the volume range 20–300 µL	5280 000 037	960001028
TS 1000, single-channel dispensing tool for the volume range 40–1,000 µL	5280 000 053	960001036
TM 10-8, 8-channel dispensing tool for the volume range 0.2–10 µL	5280 000 304	
TM 50-8, 8-channel dispensing tool for the volume range 1–50 µL	5280 000 215	960001044
TM 300-8, 8-channel dispensing tool for the volume range 20–300 µL	5280 000 231	960001052
TM 1000-8, 8-channel dispensing tool for the volume range 40–1,000 µL	5280 000 258	960001061
Holder for 6 dispensing tools	5075 774 003	960001109
epT.I.P.S.® Motion pipette tips		
Pipette tips in individual racks, for use with the epMotion®. Tip type and size are automatically recognized on the device. 96	epT.I.P.S.®/rack	
Without filter		
10 µL, sterile, volume range 0.2–10 µL, 10 × 96 tips in racks	0030 015 185	
50 µL, sterile, volume range 1–50 µL, 10 × 96 tips in racks	0030 015 207	
300 µL, sterile, volume range 20–300 µL, 10 × 96 tips in racks	0030 015 223	
1,000 µL, sterile, volume range 40–1,000 µL, 10 × 96 tips in racks	0030 015 240	
10 µL, Eppendorf Quality, volume range 0.2–10 µL, 10 × 96 tips in racks	0030 014 383	
50 µL, Eppendorf Quality, volume range 1–50 µL, 10 × 96 tips in racks	0030 014 405	
300 µL, Eppendorf Quality, volume range 20–300 µL, 10 × 96 tips in racks	0030 014 448	
1,000 µL, Eppendorf Quality, volume range 40–1,000 µL, 10 × 96 tips in racks	0030 014 480	
10 µL, Reloads, Eppendorf Quality, volume range 0.2–10 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 545	
50 µL, Reloads, Eppendorf Quality, volume range 1–50 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 421	
300 µL, Reloads, Eppendorf Quality, volume range 20–300 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 464	
1,000 µL, Reloads, Eppendorf Quality, volume range 40–1,000 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 502	
50 µL, SafeRacks, Eppendorf Quality, volume range 1–50 µL, 10 × 96 tips in racks	0030 014 600	
300 µL SafeRack, Eppendorf Quality, volume range 20–300 µL, 10 × 96 tips in racks.	0030 014 626	
1,000 µL, SafeRacks, Eppendorf Quality, volume range 40–1,000 µL, 10 × 96 tips in racks	0030 014 642	
With filter		
10 µL, sterile, PCR clean, volume range 0.2–10 µL, 10 × 96 tips in racks	0030 015 193	
50 µL, sterile, PCR clean, volume range 1–50 µL, 10 × 96 tips in racks	0030 015 215	
300 µL, sterile, PCR clean, volume range 20–300 µL, 10 × 96 tips in racks	0030 015 231	
1,000 µL, sterile, PCR clean, volume range 40–1,000 µL, 10 × 96 tips in racks	0030 015 258	
10 µL, PCR clean, volume range 0.2–10 µL, 10 × 96 tips in racks	0030 014 391	
50 µL, PCR clean, volume range 1–50 µL, 10 × 96 tips in racks	0030 014 413	
300 µL, PCR clean, volume range 20–300 µL, 10 × 96 tips in racks	0030 014 456	
1,000 µL, PCR clean, volume range 40–1,000 µL, 10 × 96 tips in racks	0030 014 499	
10 µL, Reloads, PCR clean, sterile, volume range 0.2–10 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 561	
50 µL, Reloads, PCR clean, sterile, volume range 1–50 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 529	
300 µL, Reloads, PCR clean, sterile, volume range 20–300 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 537	
1000 µL Reloads, PCR clean, sterile, volume range 40–1,000 µL, 10 × 96 tips in racks.	0030 014 570	
10 µL, Reloads, PCR clean, volume range 0.2–10 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 553	
50 µL, Reloads, PCR clean, volume range 1–50 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 430	
300 µL, Reloads, PCR clean, volume range 20–300 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 472	
1,000 µL, Reloads, PCR clean, volume range 40–1,000 µL, Reloads, 12 × 2 trays with 96 tips	0030 014 510	
50 µL, SafeRacks, PCR clean, volume range 1–50 µL, 10 × 96 tips in racks	0030 014 618	
300er SafeRacks, PCR clean, volume range 20–300 µL, 10 × 96 tips in racks	0030 014 634	
1,000 µL, SafeRacks, PCR clean, volume range 40–1,000 µL, 10 × 96 tips in racks	0030 014 650	
TipHolder, for epT.I.P.S.® Motion Reloads	5075 751 399	
Box for epT.I.P.S.® Motion 10/50/300 µL Reload, set of 10 empty boxes, reusable	0030 014 669	
Box for epT.I.P.S.® Motion 1000 µL Reload, set of 10 empty boxes, reusable	0030 014 677	
Clip for epT.I.P.S.® Motion Reloads, set of 5 clips, mandatory use to clamp trays on the Box for epT.I.P.S.® Motion 10/50/300/1000 µL Reloads, autoclavable	5075 751 070	

Description	International Ordering no.	North America Ordering no.
Other consumables for epMotion®		
10 mL reservoir , for use with reservoir rack, 5 × 10 large volume reservoirs, PCR clean	0030 126 521	
30 mL reservoir , for use with reservoir rack, 5 × 10 large volume reservoirs, PCR clean	0030 126 505	960051009
100 mL reservoir , for use with reservoir rack, 5 × 10 large volume reservoirs, PCR clean	0030 126 513	960051017
Reservoir 400 mL , large reservoir in SLAS/ANSI format gippable, 10 pieces, e.g. for use in vacuum station	5075 751 364	960002229
Low dead volume reservoir 195 mL , single well with 96 cavities, PCR clean, 4 × 5 reservoirs	0030 126 556	
Waste bags , for waste bag holder of epMotion® 5070/5073/5075 with SN < 10,000 and epMotion® 5073 with SN > 10,000. Up to 7 L capacity, ideal for disposing biologically hazardous waste, autoclavable, PP material, thickness 50 µm, transparent, 50 bags.	5075 751 763	
Waste bags , for waste bag holder of epMotion® 5070/5073/5075 with SN < 10,000 and epMotion® 5073 with SN > 10,000. Up to 7 L capacity, ideal for disposing non hazardous waste, autoclavable, PP material, thickness 50 µm, transparent, 50 bags.	5075 751 780	
Waste bags , for waste bag holder of epMotion® 5075 with SN > 10,000. Up to 10 L capacity, ideal for disposing biologically hazardous waste, autoclavable, PP material, thickness 50 µm, transparent, 50 bags.	5075 752 034	
Waste bags , for waste bag holder of epMotion® 5075 with SN > 10,000. Up to 10 L capacity, ideal for disposing non hazardous waste, autoclavable, PP material, thickness 50 µm, transparent, 50 bags.	5075 752 050	
Reservoir rack modules TC Reservoir rack modules are inserted in a reservoir rack. They can be temperature controlled with a thermal module for heating and cooling.		
PCR 0.2 mL , for 8 × 0.2 mL PCR tubes	5075 799 049	960002601
Safe-Lock Tubes , for 4 × 0.5/1.5/2 mL tubes	5075 799 081	960002620
Tubes Ø 12 mm , for 4 × Ø 12 mm tubes	5075 799 103	960002630
Tubes Ø 16 mm , for 4 × Ø 16 mm tubes	5075 799 120	960002640
Eppendorf Tubes® 5.0 mL, for 4 × 5 mL tubes	5075 799 340	
15 mL Tubes »conical« for 4 × Ø 17 mm tubes	5075 799 162	960002650
50 mL Tubes »conical« for 2 × Ø 29 mm tubes	5075 799 189	960002660
Reservoir adapter 10 mL , for 1 × epMotion® reservoir 10 mL	5075 799 421	
Reservoir adapter 30 mL , for 1 × epMotion® reservoir 30 mL	5075 799 146	960002670
Reservoir adapter 100 mL , for 1 × epMotion® reservoir 100 mL	5075 799 260	960002680
Racks For individual glass or plastic tubes, no temperature control		
Ø 12 mm × 60 mm max. length	5075 792 109	960002351
Ø 12 mm × 100 mm max. length	5075 763 001	960002059
Ø 13 mm × 60 mm max. length	5075 792 087	960002342
Ø 13 mm × 100 mm max. length	5075 762 005	960002041
Ø 14 mm × 60 mm max. length	5075 792 060	960002334
Ø 14 mm × 100 mm max. length	5075 792 001	960002369
Ø 15 mm × 60 mm max. length	5075 792 044	960002326
Ø 15 mm × 100 mm max. length	5075 792 028	960002377
Ø 16 mm × 60 mm max. length	5075 776 006	960002164
Ø 16 mm × 100 mm max. length	5075 760 002	960002032
Ø 17 mm × 60 mm max. length	5075 775 000	960002156
Ø 17 mm × 100 mm max. length	5075 761 009	960002024
Rack for 24 x HPLC tubes , Ø 12 mm × 40 mm max. length	5075 792 125	960002380
Rack for 96 x 1.5/2.0 mL screw-cap tubes , requires 2 positions on the deck	5075 791 005	960002318
Rack for 24 Safe-Lock tubes 0.5 mL/1.5 mL/2.0 mL, cannot be tempered, with adapter sleeves for 0.5 mL	5075 751 453	
Rack for 24 Safe-Lock tubes 1.5 mL/2.0 mL, no temperature control	5075 751 275	
Rack ILMN tubes , for single tubes, 40 x ø 8,4 mm and 12 x ø 11,2 mm	5075 751 747	
ReservoirRack , for ReservoirRack Modules TC and max. 7 of 10 mL, 30 mL and 100 mL reagent reservoiris	5075 754 002	960002148
ReservoirRack 3 , for max. three 10 mL, 30 mL and 100 mL reservoirs, only for 5075t	5075 754 070	
Adapter for Covaris Rack 12 x 8 Strips , for E-series rack 12 place 8 microTUBEStrip V2 from Covaris PN 500444	5075 751 801	
Thermorack Rotor/Tubes , for temperature controlled pipetting of Qiagen® Rotor-Disc® 72/100 and 20 tubes 1.5/2.0 mL	5075 751 526	
Accessories for PCR/real-time PCR and NGS applications		
TipHolder 73 enables positioning epT.I.P.S. Motion Reloadtrays above a rack of epT.I.P.S. Motion and thus doubles the tip capacity in an SLAS/ANSI position reusable and autoclavable, epBlue 40.7 or higher required	5075 751 879	
Reservoir Rack Module NGS for use with epMotion ReservoirRacks, for distributing samples and reagents for NGS-specific applications, space for 20 ILMN tubes, 4 Eppendorf Tubes 5.0 mL, 4 Eppendorf Tubes 1.5 mL	5075 751 917	
Reservoir Rack Module PCR for use with epMotion ReservoirRacks, for distributing samples and reagents for PCR, space for 48 PCR tubes 0.2 mL incl. 2 strip adapters 0.2 mL, (alternatively, space sufficient for 16 index tubes and up to 32 PCR tubes 0.2 mL), 12 Eppendorf Tubes 1.5/2.0 mL, 1 Eppendorf Tube 5.0 mL	5075 751 933	
Rack Module Tips for use with epMotion ReservoirRacks, for providing up to 16 epT.I.P.S.Motion per module and up to 7 modules per ReservoirRack (112 tips in total), set consisting of 7 modules	5075 751 950	
Tip Tool	5285 000 000	

Description	International Ordering no.	North America Ordering no.
Accessories for PCR/real-time PCR and NGS applications		
Thermoblock PCR 96 OC , for temperature control and safe piercing of sealed, semi-skirted PCR plates with orientation control	5075 751 666	
Thermoadapter for Microplate 96/V/U , for holding and temperature control of Eppendorf Microplates with V or U bottom	5075 751 577	
Rack Smart , for holding one tube rack of SmartCycler reaction tubes	5075 790 009	960002520
Rack LC , for holding up to 96 × 20 µL or 100 µL LightCycler capillaries, for use with, e.g., Centrifuge 5804/5804 R or 5810/5810 R, set of 2	5075 795 000	960002511
Thermorack CB 100 µL , for holding up to 384 x 0.1 mL strip tubes in the Corbett Research Rotor-Gene 3000	5075 767 031	960002500
Thermoadapter for PCR plates, 96-well, skirted , for use with 96 × 0.2 mL tubes or a 96-well PCR plate	5075 787 008	960002199
Thermoadapter for PCR plates, 384-well, skirted , for use with a 384-well PCR plate	5075 788 004	960002202
Thermoadapter Frosty , combination of height adapter and PCR-Cooler for cooling of skirted PCR plates	5075 789 000	960002300
Thermoblock for PCR plates, 96-well , for heating or cooling of PCR plates; plates are exchangeable using the gripper	5075 766 000	960002083
Thermoblock for PCR plates, 384-well , for heating or cooling of PCR plates; plates are exchangeable using the gripper	5075 767 007	960002091
Thermoadapter LC , designed to hold the sample cartridge of the MagNa Pure® LC System	5075 751 305	5075751305
Thermorack for 24 x 0.5 mL Safe-Lock Tubes , for a supply of 24 test tubes, temperature control	5075 769 000	960002067
Thermorack for 24 x 1.5/2 mL Safe-Lock Tubes , for a supply of 24 test tubes, temperature control	5075 771 004	960002075
Thermorack, for 24 Cryo tubes , temperature control	5075 777 055	960002491
Thermorack TMX for 24 Safe-Lock tubes, for a supply of 0.5 mL tubes on Eppendorf ThermoMixer® of 5073t, 5075t, 5075vt, temperature control no separation	5075 751 160	960002070
Thermorack TMX for 24 Safe-Lock tubes, for a supply of 1.5/2 mL tubes on Eppendorf ThermoMixer® of 5073t, 5075t, 5075vt, temperature control no separation	5075 751 186	960002080
Thermoadapter for Deep Well Plates, 96 wells/1,000 µL	5075 751 054	960002391
Thermoblock DWP 2000 for 2 mL Eppendorf Deepwell Plates 96/2000	5075 751 330	
Eppendorf Magnum FLX® Magnet Adapter , universal ring magnet plate for rapid magnetic bead separation. Supports most 96-well plates (standard, PCR, deepwell)	5075 751 836	
Accessories for vacuum station		
Vac frame holder	5075 778 009	960002237
Vac lid	5075 779 005	960002245
Mat , for vacuum lid	5075 793 008	960002407
Vac frame 1 , for skirted filter plates	5075 784 009	960002253
Vac frame 2 , for semi-skirted filter plates	5075 785 005	960002261
Collection Plate Adapter , for handling collection microtubes in microtube racks	5075 785 030	960002531
Collection Plate Adapter for MN Kits	5075 785 064	960002571
Channeling Plate Adapter 96 , for vacuum-processing foaming solutions through multiwell plates, set of 10	5075 794 004	960002540
Vac Thermo Lid , for effective drying of filter plates in epMotion® 5075 VAC vacuum chamber Note: additional thermal module required	5075 796 007	960002551
VAC Lid Sealing Mats, 5 pcs.	5075751038	960002451
Other accessories		
Gripper , incl. holder, for transporting plates on the worktable and automatic operation of the vacuum manifold	5282 000 018	960002270
Holder for Gripper , compatible with 5073 (SN < 6000) and all 5075	5075 759 004	960002211
Gripper Tower , enables positioning epT.I.P.S.® Motion Reload trays in the parking position of the gripper and holds the gripper when it is not in use; only for use with epMotion® 5073 as of serial number 6000; Gripper 5282 000 018 needs to be ordered additionally	5075 751 895	
Liquid Waste System , consisting of splash protection, waste funnel, waste tub, tub lid. For use on epMotion® 5073 and epMotion® 5075 with SN > 10,000.	5075 799 448	
PrepRack for 24 Safe-Lock tubes 2 mL	5073 751 006	
Height adapter 85 mm for uniform levels of labware; enables faster processing of the plate	5075 751 003	960002105
Height adapter 55 mm for uniform levels of labware; enables faster processing of the plate	5075 752 000	960002113
Height adapter 40 mm for uniform levels of pipette tips; enables faster processing of the plate	5075 755 009	960002121
Adapter sleeves , 1 set = 25 pieces for reconfiguration of the thermo rack (1.5/2.0 mL) for use with 0.5 mL tubes	5075 772 000	960002172
Tipholder for epMotion® reload tips	5075 751 399	

Your local distributor: www.eppendorf.com/contact
Eppendorf SE · Barkhausenweg 1 · 22339 Hamburg · Germany
eppendorf@eppendorf.com · www.eppendorf.com

www.eppendorf.com/automation

Eppendorf SE reserves the right to modify its products and services at any time. This brochure is subject to change without notice. Although prepared to ensure accuracy, Eppendorf SE assumes no liability for errors, or for any damages resulting from the application or use of this information. Referring to this brochure alone cannot as such provide for or replace reading and respecting the current version of the operating manual."

MVS® is a registered trademark of Artel Inc., USA. MagNa pure® is a registered trademark of Roche Diagnostics GmbH. Qiagen® and Rotor-Disc® are registered trademarks of Qiagen GmbH, Germany, Eppendorf®, the Eppendorf Brand Design, epMotion®, Eppendorf Magnum FLX® Eppendorf Tubes®, epT.I.P.S.®, epBlue® and Eppendorf ThermoMixer® are registered trademarks of Eppendorf SE, Germany. All rights reserved, including graphics and images. Order No.: AA01 030 320/EN1/0.1T/0323/SSO/WEID. Copyright © 2023 by Eppendorf SE.